
happybrush®

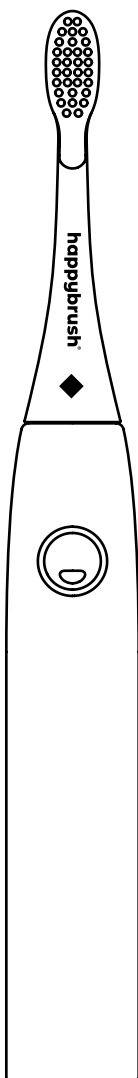ECO VIBE LITE

SCHALL-ZAHNBÜRSTE

SONIC TOOTHBRUSH | Brosse à dents sonique
CEPILLO DE DIENTES SÓNICO | SPAZZOLINO DA DENTI SONICO

DE HALLO, ICH BIN HAPPYBRUSH!

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SUPPORT

happybrush GmbH
Lindwurmstr. 76
80337 München

Hotline: 089 24412635
E-Mail: support@happybrush.de
Web: www.happybrush.de

EN HELLO, I'M HAPPYBRUSH!

HIGH-PERFORMANCE POWER BATTERY

My power batter (Li-Ion) will accompany you on your adventures for up to 6 weeks.

100% RECYCLED MATERIAL

My handle is made from 100% recycled material.

PRACTICAL “LOCK” FUNCTION

You can simply “lock” me when you travel (5 second power button).

HIGH-QUALITY TRIPLE BRISTLES

I have triple-twisted bristles which pamper your teeth with their special surface. My indicator bristles show when it is time to replace your brush head.

EFFECTIVE SONIC TECHNOLOGY

My sonic technology ensures effective and gentle cleaning with up to 25,000 vibrations.

INTEGRATED 2 MINUTE TIMER

My integrated 2-minute quadrant timer alerts you every 30 seconds to change your brushing zone.

SAFETY INSTRUCTIONS

EXPLANATION OF SYMBOLS

Safety instructions: Read through these carefully and follow them to prevent personal injury and damage to property.

For indoor use only

USB output

Read operating instructions

Dispose of the device according to the WEEE directive.

- This device may only be used for its intended purpose of tooth cleaning and oral care. Commercial use is excluded.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Children shall not play with the appliance.
- Cleaning and user maintenance shall not be made by children without supervision
- If the electric toothbrush or the mains cable of the device is damaged, it must be replaced by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid

any hazards. This appliance contains batteries that are non-replaceable.

- Only use the toothbrush with the provided mains cable.
- **WARNING:** Only use a certified safety extra-low voltage (SELV) mains plug with an output voltage of 5V and 1A to charge the toothbrush. For safe charging in a damp environment (e.g. in the bathroom), only use a splash-proof mains plug with protection class IPX4. The use of a non-certified mains plug can lead to danger or serious injury.

Danger of Electric Shock from Water!

- Never immerse the mains cable in water or other liquids and ensure that it does not fall into water or become wet.
- If the toothbrush and mains cable fall into water during charging, immediately disconnect the power supply. Do not attempt to pull the toothbrush and mains cable out of the water while they are connected to the mains power!
- Never touch the mains cable or mains plug with wet hands when these components are connected to the mains power.

Risk of Damage to Material and Property

- Position the mains cable in such a way that it does not become squashed or bent and does not come into contact with hot surfaces.
- Remove the mains plug from the plug socket if a fault occurs during charging or before a thunderstorm.
- When you pull it out of the plug socket, always pull on the mains plug and never the mains cable!

Warning of fire!

- Do not use the device in rooms containing easily ignitable dust or poisonous and explosive fumes.
- Do not cover the charging station during charging.

The Proper Handling of the Battery

- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected area with abundant fresh water immediately and contact a doctor as soon as possible.
- The battery must not be taken apart, thrown into a fire, immersed in liquids or short-circuited.

PREPARATION AND INSTRUCTIONS

CHARGING

1. Connect the plug of the USB-C cable to a suitable mains plug (5V 1A).
2. Open the cover of the charging port at the bottom end of the handle. Connect the other end of the USB-C cable directly to your happybrush to start the charging process. A full charge of the battery takes about 10 hours. The first charging process should take 24 hours. During charging, the LED will pulse. The charging process is complete when the LED lights up continuously. It is normal for the handle to warm up during charging.

Important: After finishing the charging process, make sure to fully close the charging port cover to prevent water from entering.

BRUSHING INSTRUCTIONS

1. Moisten my bristles with a bit of water and then apply toothpaste - happybrush toothpaste works best. Now you can start brushing!
2. Place my brush head at a 45-degree angle to your teeth and guide me over your teeth with stroking movements and only light pressure on the gums.
3. Start at the end of your upper dental arch and then run me over the outside of your teeth, then the inside. Repeat the same for your lower dental arch.
4. The integrated quadrant timer gives you a short signal every 30 seconds (1x vibration) so that you know when it's time to change the brushing zones.
5. Clean your chewing surfaces with forward and backward movements. After 2 minutes, the regular brushing cycle is over. Your happybrush will then give you a slightly longer signal (2x vibration).

MODE OF OPERATION


READY, SET...

Use the power button to start me.

GO!

I vibrate every 30 seconds so that you can equally brush your 4 quadrants.. You can pause or stop me by pressing my power button.

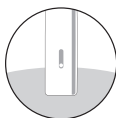ALREADY DONE

After 2 minutes, I make two brief vibrations. You can keep on brushing if you want. My red LED flashes after switching off when my battery level drops below 15 percent.

ON THE ROAD

You can read my battery level by pressing the power button for 2 seconds. And if you press the power button for 5 seconds, I'll "lock" for travel.

CLEANING

HANDLE

1. Remove the brush head and rinse the metal rod with warm water. Please do not press on the seal around the metal rod with sharp-edged objects, as this could damage it.
2. Simply wipe the handle with a dry cloth.

BRUSH HEAD

1. Rinse the brush head and bristles under running water after each use. Then, lightly shake out the water from the brush head.
2. Remove the brush head from the handle and rinse the connection between the handle and brush head with warm water.
3. For optimal results, use a new brush head every 2-3 months. The indicator bristles will show when it's time: Simply replace it when the smile fades.

Important: The handle and the brush head should not be cleaned in the dishwasher. Please avoid using harsh cleaning agents, brushes with metal or nylon bristles, or sharp or metallic cleaning tools such as knives or similar items.

DISPOSAL

RECYCLING

The environment is very important to us. When you no longer use your happybrush, do not throw it in with your normal household waste.

Please note that your happybrush contains an integrated Li-ion rechargeable battery. Dispose of your happybrush at a recycling centre for used electrical and electronic devices. More information is available from your local authorities.

To protect the environment, non-rechargeable and rechargeable batteries must not be disposed of with normal household waste, but must be taken to suitable collection points. If the device is disposed of, the integrated battery may only be removed by qualified specialists. Also note the relevant statutory regulations governing the disposal of batteries.

You can also return your happybrush directly to us after use and we will take care of its proper disposal. Please simply contact us at support@happybrush.de

NOTES ON POSSIBLE FAULTS

PROBLEM	POSSIBLE CAUSE	SOLUTION
Your happybrush does not start	The battery of your happybrush is not charged.	Charge your happybrush. The charge level may vary before the first use. Charge for approx. 24 hours before first use.
	The travel lock is activated.	Press and hold the power button for 5 seconds to release the travel lock.
The cleaning performance of your happybrush decreases	Your happybrush is not fully charged.	Charge your happybrush. The charge level may vary before the first use. Charge for approx. 24 hours before first use.
	The brush head has been used for too long.	We recommend changing the brush heads regularly (every 2-3 months).
The LEDs do not indicate the charging process.	The cable is not correctly connected.	Check that your cable is correctly connected via USB-C or to the mains plug and that the mains plug is correctly inserted into the socket.
Your happybrush switches itself off while you are brushing your teeth.	The battery may be completely discharged.	Recharge your happybrush for up to 24 hours.
Your happybrush is louder than usual.	Your brush head needs to be replaced or is defective.	Change your brush head.

PROBLEM	POSSIBLE CAUSE	SOLUTION
The brush head does not fit on the happybrush.	You are using the wrong brush head.	Look out for the name ECO VIBE 3 or 3+ when buying the brush heads
	The inner part of the old brush head is still attached to the metal pin.	Carefully detach the inner part from the metal pin of your happybrush.

NOTES ON THE PURCHASE OF BRUSH HEADS

Don't worry, we always have supplies for you! You can find the right brush head for you either at the retail partner where you previously purchased your happybrush products or in our online shop, where you will also find an overview of all sales outlets. Look out for the product name ECO VIBE 3 or 3+ when buying the brush heads. We recommend changing the brush attachments regularly (every 2-3 months) to ensure thorough and hygienic cleaning.

TECHNICAL SPECIFICATIONS

Battery	3.7V  1200 mAh (Li-Ion)
Input	5V  1A
Protection level	IPX7 (handle)
Vibrations	25.000 per minute
Type	HBV3ES
Timer	4-quadrants timer (every 30 seconds)
Operation time	approx. 6 weeks at 2 x 2 minutes per day
Charging time	Charge for 24 hours the first time, thereafter 10 hours

 The CE label on this product is our way of confirming that it complies with the valid, product-specific European guidelines.

WARRANTY

Your happybrush comes with a 1 year warranty as of the date of purchase. If defects of material or function should occur during the 1 year warranty period, we will repair or replace your happybrush.

The original proof of purchase (original receipt or payment confirmation) is to be sent to us as the warranty and proof of ownership. The receipt must show the date of purchase, your name and the full product type description of the device. The warranty applies only to the original purchaser. Should the device be sold to a second buyer, the warranty is no longer valid.

All claims from this warranty (which includes both the primary claim as well as any secondary claims after a repair or a replacement) are subject to a 12 month limitation period. The limitation period for the primary warranty starts from the moment the defect is known or gross negligent ignorance of the defect. The limitation period for secondary claims begins with the delivery of the repaired or replaced product.

The product may be replaced with the same model or a comparably equipped model. The replacement or repair of the product does not justify a new warranty. The warranty does not apply to devices or device parts subject to normal wear and tear and thus viewed as wear parts. The warranty is void if the device has been damaged, not used properly or not maintained properly. Repairs may only be done by happybrush GmbH. To ensure proper utilization of the device, all of the instructions which are listed in the operating instructions must be complied with precisely. Using it and handling it in ways that are not recommended in the instructions or for which warnings exist are to be absolutely avoided.

Important: The warranty and warranty period for specific devices and materials can differ from those shown here. These warranty restrictions have no impact on the statutory claim to warranty. The following device parts do not fall under the guarantee: Brush heads; damage due to misuse, wrong use, negligence or changing the product; normal wear and tear such as chipped spots, scratches, abrasions, discolourations or fading.

SUPPORT

happybrush GmbH
Lindwurmstr. 76
80337 Munich

Hotline: +49 (0)89 24412635
Email: support@happybrush.com
Web: www.happybrush.com

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