
# happybrush®
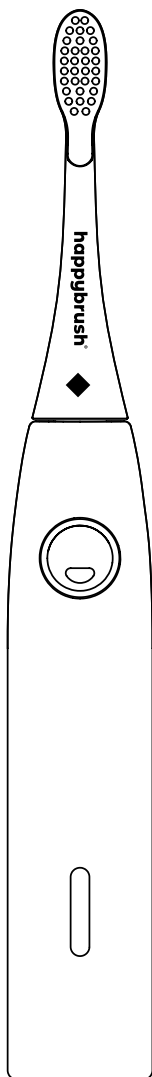
## ECO VIBE 3

### SCHALL-ZAHNBÜRSTE

SONIC TOOTHBRUSH | Brosse à dents sonique  
CEPILLO DE DIENTES SÓNICO | SPAZZOLINO DA DENTI SONICO

# DE HALLO, ICH BIN HAPPYBRUSH!

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## SUPPORT

happybrush GmbH  
Lindwurmstr. 76  
80337 München

Hotline: 089 24412635  
E-Mail: [support@happybrush.de](mailto:support@happybrush.de)  
Web: [www.happybrush.de](http://www.happybrush.de)

# EN HELLO, I'M HAPPYBRUSH!

## HIGH-PERFORMANCE POWER BATTERY

My power batter (Li-Ion) will accompany you on your adventures for up to 6 weeks.

## 100% RECYCLED MATERIAL

My handle and charging station housing are made from 100% recycled material.

## CUSTOM CLEANING MODES

NORMAL: Thorough and effective cleaning

SENSITIVE: Gentle cleaning of sensitive teeth

POLISH: Special polishing of individual teeth

## COMPACT CHARGING STATION

You can charge me in many places with my practical USB-C cable.

## HIGH-QUALITY TRIPLE BRISTLES

I have triple-twisted bristles which pamper your teeth with their special surface. My indicator bristles show when it is time to replace your brush head.

## EFFECTIVE SONIC TECHNOLOGY


My sonic technology ensures effective and gentle cleaning with up to 40,000 vibrations.

## INTEGRATED 2 MINUTE TIMER

My integrated 2-minute quadrant timer alerts you every 30 seconds to change your brushing zone.

## PRACTICAL "LOCK" FUNCTION

You can simply "lock" me when you travel (5 second power button).
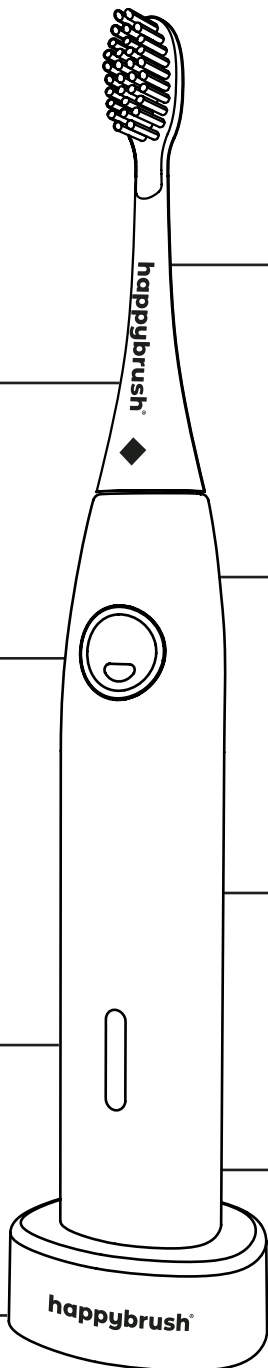
# SAFETY INSTRUCTIONS

## EXPLANATION OF SYMBOLS

Safety instructions: Read through these carefully and follow them to prevent personal injury and damage to property.

 For indoor use only

USB output

Read operating instructions

Dispose of the device according to the WEEE directive.

- This device may only be used for its intended purpose of tooth cleaning and oral care. Commercial use is excluded.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved.
- Cleaning and user maintenance shall not be made by children without supervision
- Children shall not play with the appliance.
- The charging station must always be disconnected from the mains power when left unattended and before cleaning it.
- Always check the charging station for damage before putting it into operation. Only use it if the charging station and

the mains cable have no signs of damage!

- If the electric toothbrush, the charging station or the mains cable of the device is damaged, it must be replaced by the manufacturer, customer service department or a similarly qualified person (e.g. specialist workshop) in order to avoid any hazards. This appliance contains batteries that are non-replaceable.
- Only use the toothbrush with the provided charging station.
- Warning: Only use a certified safety extra-low voltage (SELV) mains plug with an output voltage of 5V and 1A to charge the toothbrush. For safe charging in a damp environment (e.g. in the bathroom), only use a splash-proof mains plug with protection class IPX4. The use of a non-certified mains plug can lead to danger or serious injury.

### **Danger of Electric Shock from Water!**

- Never immerse the charging station or mains cable in water or other liquids and ensure that these components cannot fall into water or become wet.
- Use the charging station only in closed rooms.
- Should the charging station fall into water, cut off the power supply immediately. Never attempt to pull the charging station out of the water while it is still connected to the mains power!
- Never touch the charging station, mains cable or mains plug with wet hands when these components are connected to the mains power.

## Risk of Damage to Material and Property

- Position the mains cable in such a way that it does not become squashed or bent and does not come into contact with hot surfaces.
- Remove the mains plug from the plug socket if a fault occurs during charging or before a thunderstorm.
- When you pull it out of the plug socket, always pull on the mains plug and never the mains cable!

## Warning of fire!

- Do not use the device in rooms containing easily ignitable dust or poisonous and explosive fumes.
- Do not cover the charging station during charging.

## The Proper Handling of the Battery

- If the battery has leaked, do not allow your skin, eyes and mucous membranes to come into contact with the battery acid. If you touch battery acid, rinse the affected area with abundant fresh water immediately and contact a doctor as soon as possible.
- The battery must not be taken apart, thrown into a fire, immersed in liquids or short-circuited.

# PREPARATION AND INSTRUCTIONS

## CHARGING

1. Connect the plug of the USB-C cable to a suitable mains plug (5V 1A).
2. Place your happybrush onto the charging station. A complete charge of the battery takes approx. 10 hours. The first charge should take 24 hours. A pulsating LED bar indicates that charging is in progress. Charging is complete when all the LEDs light up permanently. It is normal for the handle to become warm during charging.

Important: Charging occurs via induction and can be affected by metallic objects near the charging station. Keep at least 10 cm between the charging station and any metallic objects.


## BRUSH HEADS

Insert the brush head onto the metal tip of the handle until it sits firmly.

## BRUSHING INSTRUCTIONS

1. Moisten my bristles with a bit of water and then apply toothpaste - happybrush toothpaste works best. Now you can start brushing!
2. Place my brush head at a 45-degree angle to your teeth and guide me over your teeth with stroking movements and only light pressure on the gums.
3. Start at the end of your upper dental arch and then run me over the outside of your teeth, then the inside. Repeat the same for your lower dental arch.
4. The integrated quadrant timer gives you a short signal every 30 seconds (1x vibration) so that you know when it's time to change the brushing zones.
5. Clean your chewing surfaces with forward and backward movements. After 2 minutes, the regular brushing cycle is over. Your happybrush will then give you a slightly longer signal (2x vibration). You can now polish areas with discolouration.

## MODE OF OPERATION
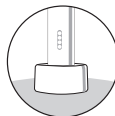
### CHARGING

My 4 LEDs pulsate during charging. If all the LEDs are lit up, it means I am fully charged.

### READY, SET...

Use the power button to start me and change operating modes simply by pressing the button again (within 5 seconds).

### GO!

I vibrate every 30 seconds so that you can equally brush your 4 quadrants.. You can pause or stop me by pressing my power button. My LEDs will then light up and down.

### ALREADY DONE

After 2 minutes, I make two brief vibrations. You can keep on brushing if you want. My bottom LED flashes after switching off when my battery level drops below 10 percent.

### ON THE ROAD

You can read my battery level by pressing the power button for 2 seconds. And if you press the power button for 5 seconds, I'll "lock" for travel.

# CLEANING

## HANDLE

1. Take the brush head off and rinse the metal tip with warm water. Please do not try to clean around the seal of the metal tip with a sharp object because it could damage the seal.
2. Now, simply wipe the handle with a dry cloth.

## BRUSH HEAD

1. Rinse the brush head and bristles under running water after every use. Then just shake excess water off the brush head.
2. Remove the brush head from the toothbrush and clean the connection between the handle and the brush head with warm water.
3. Use a new brush head every 2-3 months for best results. The indicator bristles will tell you when: Simply change it when the smile fades.

## CHARGING STATION

1. Remove the plug from the mains socket before cleaning the charging station.
2. Then wipe the charging station with a slightly damp cloth.

Important: The handle, the brush heads, and the charging station are not to be cleaned in the dishwasher. Do not submerge the charging station in water. Do not use abrasive cleaning agents. Also do not try to clean the charging station with any kind of sharp metallic object or nylon bristles.

# DISPOSAL

## RECYCLING

The environment is very important to us. When you no longer use your happybrush, do not throw it in with your normal household waste.

Please note that your happybrush contains an integrated Li-ion rechargeable battery. Dispose of your happybrush at a recycling centre for used electrical and electronic devices. More information is available from your local authorities.

To protect the environment, non-rechargeable and rechargeable batteries must not be disposed of with normal household waste, but must be taken to suitable collection points. If the device is disposed of, the integrated battery may only be removed by qualified specialists. Also note the relevant statutory regulations governing the disposal of batteries.

You can also return your happybrush directly to us after use and we will take care of its proper disposal. Please simply contact us at [support@happybrush.de](mailto:support@happybrush.de)

# NOTES ON POSSIBLE FAULTS

PROBLEM	POSSIBLE CAUSE	SOLUTION
Your happybrush does not start	The battery of your happybrush is not charged.	Charge your happybrush. The charge level may vary before the first use. Charge for approx. 24 hours before first use.
	The travel lock is activated.	Press and hold the power button for 5 seconds to release the travel lock or place your happybrush on the charging station.
The cleaning performance of your happybrush decreases	Your happybrush is not fully charged.	Charge your happybrush. The charge level may vary before the first use. Charge for approx. 24 hours before first use.
	The brush head has been used for too long.	We recommend changing the brush heads regularly (every 2-3 months).
The LEDs do not indicate the charging process.	Charging station is not connected.	Check that your cable is correctly connected via USB-C or to the mains plug and that the mains plug is correctly inserted into the socket.
	Your happybrush is not placed correctly on the charging station.	Place your happybrush in the centre of the charging station.
Your happybrush switches itself off while you are brushing your teeth.	The battery may be completely discharged.	Recharge your happybrush for up to 24 hours.
Your happybrush is louder than usual.	Your brush head needs to be replaced or is defective.	Change your brush head.

PROBLEM	POSSIBLE CAUSE	SOLUTION
The brush head does not fit on the happybrush.	You are using the wrong brush head.	Look out for the name ECO VIBE 3 or 3+ when buying the brush heads
	The inner part of the old brush head is still attached to the metal pin.	Carefully detach the inner part from the metal pin of your happybrush.

## NOTES ON THE PURCHASE OF BRUSH HEADS

Don't worry, we always have supplies for you! You can find the right brush head for you either at the retail partner where you previously purchased your happybrush products or in our online shop, where you will also find an overview of all sales outlets. Look out for the product name ECO VIBE 3 or 3+ when buying the brush heads. We recommend changing the brush attachments regularly (every 2-3 months) to ensure thorough and hygienic cleaning.

## TECHNICAL SPECIFICATIONS

<b>Battery</b>	3.7V  1200 mAh (Li-Ion)
<b>Input</b>	5V  1A
<b>Protection level</b>	IPX7
<b>Vibrations</b>	Normal: 25.000 per minute Sensitive: 30.000 per minute Polishing: 30.000-40.000 per minute
<b>Type</b>	HBV3
<b>Timer</b>	4-quadrants timer (every 30 seconds)
<b>Operation time</b>	approx. 6 weeks at 2 x 2 minutes per day
<b>Charging time</b>	Charge for 24 hours the first time, thereafter 10 hours

The CE label on this product is our way of confirming that it complies with the valid, product-specific European guidelines.

## WARRANTY

Your happybrush comes with a 1 year warranty as of the date of purchase. If defects of material or function should occur during the 1 year warranty period, we will repair or replace your happybrush.

The original proof of purchase (original receipt or payment confirmation) is to be sent to us as the warranty and proof of ownership. The receipt must show the date of purchase, your name and the full product type description of the device. The warranty applies only to the original purchaser. Should the device be sold to a second buyer, the warranty is no longer valid.

All claims from this warranty (which includes both the primary claim as well as any secondary claims after a repair or a replacement) are subject to a 12 month limitation period. The limitation period for the primary warranty starts from the moment the defect is known or gross negligent ignorance of the defect. The limitation period for secondary claims begins with the delivery of the repaired or replaced product.

The product may be replaced with the same model or a comparably equipped model. The replacement or repair of the product does not justify a new warranty. The warranty does not apply to devices or device parts subject to normal wear and tear and thus viewed as wear parts. The warranty is void if the device has been damaged, not used properly or not maintained properly. Repairs may only be done by happybrush GmbH. To ensure proper utilization of the device, all of the instructions which are listed in the operating instructions must be complied with precisely. Using it and handling it in ways that are not recommended in the instructions or for which warnings exist are to be absolutely avoided.

Important: The warranty and warranty period for specific devices and materials can differ from those shown here. These warranty restrictions have no impact on the statutory claim to warranty. The following device parts do not fall under the guarantee: Brush heads; damage due to misuse, wrong use, negligence or changing the product; normal wear and tear such as chipped spots, scratches, abrasions, discolourations or fading.

## SUPPORT

happybrush GmbH  
Lindwurmstr. 76  
80337 Munich

Hotline: +49 (0)89 24412635  
Email: [support@happybrush.com](mailto:support@happybrush.com)  
Web: [www.happybrush.com](http://www.happybrush.com)