beurer

FC 41

DE	Porenreiniger Gebrauchsanweisung	2	IT	Dispositivo per la pulizia dei pori Istruzioni per l'uso	56
EN	Power deep pore cleanser Instructions for use.....	16	TR	Gözenek temizleyici Kullanım kılavuzu.....	70
FR	Nettoyeur de pores Mode d'emploi.....	29	RU	Очиститель пор Инструкция по применению	83
ES	Purificador de poros Instrucciones de uso.....	42	PL	Oczyszczacz do porów Instrukcja obsługi	96

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Read these instructions for use carefully and keep them for later use, be sure to make them accessible to other users and observe the information they contain.

WARNING

- The device is only intended for domestic/private use, not for commercial use.
- This device is not suitable for persons under the age of 18.
- This device may be used by people with reduced physical, sensory or mental skills or a lack of experience or knowledge, provided that they are supervised or have been instructed on how to use the device safely and are fully aware of the consequent risks of use.
- Children must not play with the device.
- Cleaning and user maintenance must not be performed by children unless supervised.
- If the mains connection cable of this device is damaged, it must be disposed of. If it cannot be removed, the device must be disposed of.
- Under no circumstances should you open or repair the device yourself, as faultless functionality can no longer be guaranteed thereafter. Failure to adhere to this instruction will void the warranty.

Dear customer,

Thank you for choosing a product from our range. Our name stands for high-quality, thoroughly tested products for applications in the areas of heat, gentle therapy, blood pressure/diagnosis, weight, massage and air.

With kind regards,
Your Beurer team

Table of contents

1. Included in delivery	17	8.1 Carrying out a skin test	23
2. Signs and symbols	18	8.2 Charging the battery	24
3. Getting to know your device	19	8.3 Please note the following before use.....	24
4. Intended use.....	19	8.4 Application areas.....	25
5. Warnings and safety notes	19	9. Cleaning and maintenance.....	26
6. Device description.....	21	10. Replacement parts and wearing parts	26
6.1 Device	21	11. Disposal.....	26
6.2 Attachments	22	12. Technical specifications	27
7. Initial use	22	13. Warranty/service.....	27
8. Use	23		

1. Included in delivery

Check that the exterior of the delivery cardboard packaging is intact and make sure that all components are present. Before use, ensure that there is no visible damage to the device or accessories and that all packaging material has been removed. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

- 1 x FC 41
- 1 x round attachment (large)
- 1 x round attachment (small)
- 1 x oval attachment
- 10 x replacement filters
- 2 x sealing rings
- 1 x filter cover
- 1 x USB charging cable
- 1 x set of instructions for use

2. Signs and symbols

The following symbols appear in these instructions for use and on the device:

	DANGER This device must not be used near to water or in water (e.g. in a wash basin, or in the shower or bath) – risk of electric shock!
	WARNING Warning notice indicating a risk of injury or damage to health
	IMPORTANT Safety note indicating possible damage to the device/accessory
	Note Note on important information
	Read the instructions
	Direct current The device is suitable for use with direct current only
	Disposal in accordance with the Waste Electrical and Electronic Equipment EC Directive – WEEE
	Manufacturer
	Only use indoors
	Dispose of packaging in an environmentally friendly manner
	This product satisfies the requirements of the applicable European and national directives.
	The products demonstrably meet the requirements of the Technical Regulations of the EAEU.

3. Getting to know your device

The FC 41 power deep pore cleanser enables deep-pore cleansing thanks to state-of-the-art vacuum technology. Using the round attachments, you can effectively remove spots, blackheads and excess skin flakes. Furthermore, the power deep pore cleanser has an oval attachment that you can use to improve skin elasticity and the circulation in the face.

The power deep pore cleanser has a clear LCD display for displaying the battery state and intensity level. With a total of 5 intensity levels available, the power deep pore cleanser enables personalised cleansing and care according to your needs.

4. Intended use

The device is only to be used on the face (apart from the eyelids and lips) of human adults. Do not use the device on animals!

The device is only intended for private use, and not for medical or commercial purposes. Use the device for the purpose it was developed for only and in the manner specified in these instructions for use. Any form of improper use can be dangerous. The manufacturer is not liable for damage resulting from improper or incorrect use. Do not use the device if one or several of the following warnings and safety notes applies to you. If you are unsure whether the device is suitable for you, please consult your doctor.

5. Warnings and safety notes

WARNING

- Danger of suffocation by packaging material. Keep packaging material away from children.
- Before use, inspect the device and accessories for visible damage. If you have any doubts, do not use the device and contact your retailer or the specified Customer Service address.

WARNING

Do **NOT** use the device

- on cracked skin
- on open wounds

- in the eye area or on the mouth
- on any part of the body other than the face (apart from eye area and mouth)
- during pregnancy
- if you have had facial cosmetic surgery
- if you suffer from any type of heart disease
- if you have a skin disease (e.g. herpes, vascular diseases or skin growths)
- if you have skin cancer
- if you suffer from skin irritation
- immediately after sunbathing or if you have sunburn
- if you are taking medication that contains steroids
- if you are taking blood thinners
- if you have multiple sclerosis, diabetes or an autoimmune disorder

IMPORTANT

Only specialist personnel may perform repairs on electrical devices. Improper repairs may subject users to considerable danger. In the event of a fault or damage, have the device repaired by a certified specialist workshop.

WARNING

- Ensure that the USB charging cable and mains adapter do not come into contact with water, steam or other liquids.
- Only use the device in dry rooms (e.g. never in the bath, sauna).
- Do not drop the device or subject it to any impacts.
- Only use the device with dry hands.


Notes on handling batteries

- If your skin or eyes come into contact with fluid from a rechargeable battery cell, flush out the affected area with water and seek medical assistance.
- Choking hazard! Small children may swallow and choke on batteries. Store batteries out of the reach of small children.
- Protect batteries from excessive heat.
- Risk of explosion! Never throw batteries into a fire.
- Do not disassemble, split or crush the batteries.
- Only use chargers specified in the instructions for use.
- Batteries must be charged correctly prior to use. The instructions from the manufacturer and the specifications in these instructions for use regarding correct charging must be observed at all times.

- Fully charge the battery prior to initial use.
- In order to achieve as long a battery service life as possible, fully charge the battery at least twice per year.

6. Device description

6.1 Device

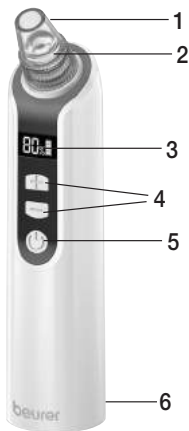1. Attachment

2. Suction head with filter cover

3. Intensity and battery state indicator

4. +/- button

For adjusting the intensity level (20%, 40%, 60%, 80% or 99%)

5. ON/OFF button

6. Micro USB connection


For charging the battery with the micro USB charging cable

6.2 Attachments


	<p>Round attachment (large) Suitable for pore cleansing on large areas. Can be used anywhere on the face, apart from on the skin around the eyes.</p>
	<p>Round attachment (small) Suitable for cleansing the pores on hard-to-reach areas (e.g. on the nose). Can be used anywhere on the face, even on sensitive skin around the eyes.</p>
	<p>Oval attachment Promotes skin elasticity and circulation in the face. Particularly suitable for use on the outside of the nostrils and around the eyes and mouth. We recommend that you use the oval attachment with intensity level 1.</p>

7. Initial use


1. Turn the attachment to take it off the device and remove the cover from the suction head.

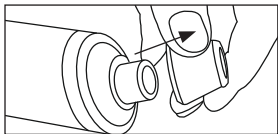
2. Place a filter on the suction head and then place the cover back over it. Make sure that the filter remains dry.

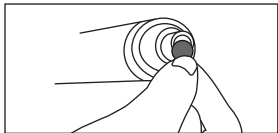
3. Place the attachment you wish to use on the suction head.

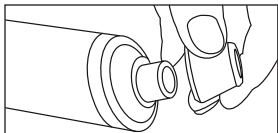4. Press the ON/OFF button to now switch on the device. The intensity and battery state indicator begins to light up.

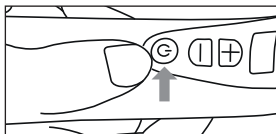5. Select your desired intensity level using the -/+ buttons. We recommend you start using a low intensity level during the first few applications so that you get used to the device gradually.

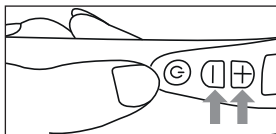
6. Move the attachment over your face. Make sure that you do not move the attachment over one part of the face more than once. Further details on treatment can be found in Section "8.4 Application areas".


i Note

The pressure you apply when using the device should be comfortable and not painful. Some slight reddening should be expected after use.

7. To end use, press the ON/OFF button. The intensity and battery state indicator goes out.

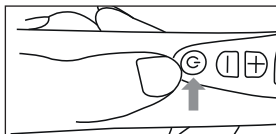8. Change the filter after each use. To do so, turn the attachment to lift it off; then remove the cover below it. Use tweezers to remove the filter.

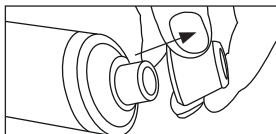8. Use

8.1 Carrying out a skin test

Carry out a skin test with the device 24 hours before use. Test the device at the lowest intensity level on small, inconspicuous areas of your neck. If, after 24 hours, the tested area is red or painful, do not use the device.

IMPORTANT

Do not use the device on the same area of skin for a long period of time in order to prevent bruises.

8.2 Charging the battery

Before using the device for the first time, let it charge for a minimum of 5 hours first. Proceed as follows:

- Connect the USB charging cable with a suitable mains adapter (output: 5V; 0.5A) and the device. Then insert the mains adapter into a suitable socket. Alternatively, you can also charge the device with the USB charging cable via the USB port on your computer or laptop.

The battery state indicator in the display displays the current battery state during charging. Once the device is fully charged,  appears in the display. The battery will last for approx. 90 minutes of operation (at the highest intensity level) when fully charged. If the battery charge indicator  flashes in the display, you have to recharge the device.

8.3 Please note the following before use

Note

Before you start using the device, make sure that your face is clean, dry and free of any traces of products such as powder, cream or make-up.

WARNING

- Do not use the device more than once a week.
- Immediately after use, the treated skin on your face may be slightly reddened. This is normal and the reddening will disappear by itself after 30 minutes at the latest.
- For one day before and one day after use, avoid activities that cause you to perspire, such as saunas, tanning beds and sun exposure.
- After use: Clean your face with water and protect it against sun exposure.

Note

The device can also be used twice a week for non-sensitive problem skin (with a break of at least 48 hours between use).

8.4 Application areas


Ensure that the part of your face you are treating is taut when you move the device over the skin. You can do this by pulling the skin taut to the side or downwards with your free hand.

Do not use the device on the same area of skin for a long period of time.


Forehead:

Move the device slowly and evenly outwards from the centre of the forehead. Make sure that you do not treat any part of your skin more than once.

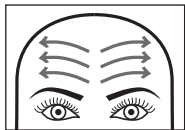Temples:

Move the device outwards towards your hairline. Repeat this on the left and right of the nose. Ensure that the device is at least 1 cm away from your eyes.

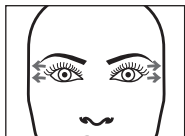Nose:

Move the device down the centre of your nose starting from the top of the bridge. Repeat this on the left and right of the nose.


Cheeks:

Move the device outwards across your cheeks.

Around the mouth:

Move the device along your upper lip from the centre to the left and right. Move the device upwards from your bottom lip to the left and right. Make sure that you do not touch your lips.

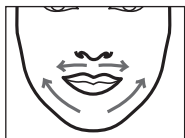9. Cleaning and maintenance

- Never submerge the device in water or other liquids.
- Disinfect the attachments with a standard disinfectant after each use.
- Change the filter after each use.
- Use a dry cloth to clean the device of any dirt. Do not use any abrasive cleaning products.
- Protect the device from impacts, moisture, dust, chemicals, marked temperature fluctuations and nearby sources of heat (ovens, heaters).

10. Replacement parts and wearing parts

Replacement parts and wearing parts are available from the corresponding service address (according to the service address list).

Please state the corresponding order number.

Replacement part	Item number and/or order number
10 replacement filters	164083
2 sealing rings	164084
Round attachment (large)	164080
Round attachment (small)	164081
Oval attachment	164082
USB charging cable	164079
Filter cover	164085

11. Disposal

For environmental reasons, do not dispose of the device in the household waste at the end of its useful life. Dispose of the device at a suitable local collection or recycling point in your country.

Dispose of the device in accordance with EC Directive – WEEE (Waste Electrical and Electronic Equipment). If you have any questions, please contact the local authorities responsible for waste disposal.

12. Technical specifications

Supply	5 V --- 0.5 A
Dimensions	approx. 169 x 43 x 43 mm
<u>Battery:</u>	
Capacity	900 mAh
Nominal voltage	3.7 V
Type designation	Lithium-ion
Weight	~150 g

13. Warranty/service

Beurer GmbH, Söflinger Straße 218, 89077 Ulm, Germany (hereinafter referred to as “Beurer”) provides a warranty for this product, subject to the requirements below and to the extent described as follows.

The warranty conditions below shall not affect the seller’s statutory warranty obligations which ensue from the sales agreement with the buyer.

The warranty shall apply without prejudice to any mandatory statutory provisions on liability.

Beurer guarantees the perfect functionality and completeness of this product.

The worldwide warranty period is 3 years, commencing from the purchase of the new, unused product from the seller.

The warranty only applies to products purchased by the buyer as a consumer and used exclusively for personal purposes in the context of domestic use.

German law shall apply.

During the warranty period, should this product prove to be incomplete or defective in functionality in accordance with the following provisions, Beurer shall carry out a repair or a replacement delivery free of charge, in accordance with these warranty conditions.

If the buyer wishes to make a warranty claim, they should approach their local retailer in the first instance: see the attached “International Service” list of service addresses.

The buyer will then receive further information about the processing of the warranty claim, e.g. where they can send the product and what documentation is required.

A warranty claim shall only be considered if the buyer can provide Beurer, or an authorised Beurer partner, with

- a copy of the invoice/purchase receipt, and
- the original product.

The following are explicitly excluded from this warranty:

- deterioration due to normal use or consumption of the product;
- accessories supplied with this product which are worn out or used up through proper use (e.g. batteries, rechargeable batteries, cuffs, seals, electrodes, light sources, attachments and nebuliser accessories);
- products that are used, cleaned, stored or maintained improperly and/or contrary to the provisions of the instructions for use, as well as products that have been opened, repaired or modified by the buyer or by a service centre not authorised by Beurer;
- damage that arises during transport between manufacturer and customer, or between service centre and customer;
- products purchased as seconds or as used goods;
- consequential damage arising from a fault in this product (however, in this case, claims may exist arising from product liability or other compulsory statutory liability provisions).

Repairs or an exchange in full do not extend the warranty period under any circumstances.

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